


Give every newborn the best start in life



BNG MEDICAL INSTRUMENTS CO.,LTD.

📍 3rd Floor, Block 3, No.14, Zhongxing Road, Xiuxin Community,
Kengzi Street, Pingshan District, Shenzhen, China
✉ Email: info@bngmed.com ☎ Tel: +86-0755-27392867
🌐 Web: www.bngmed.com

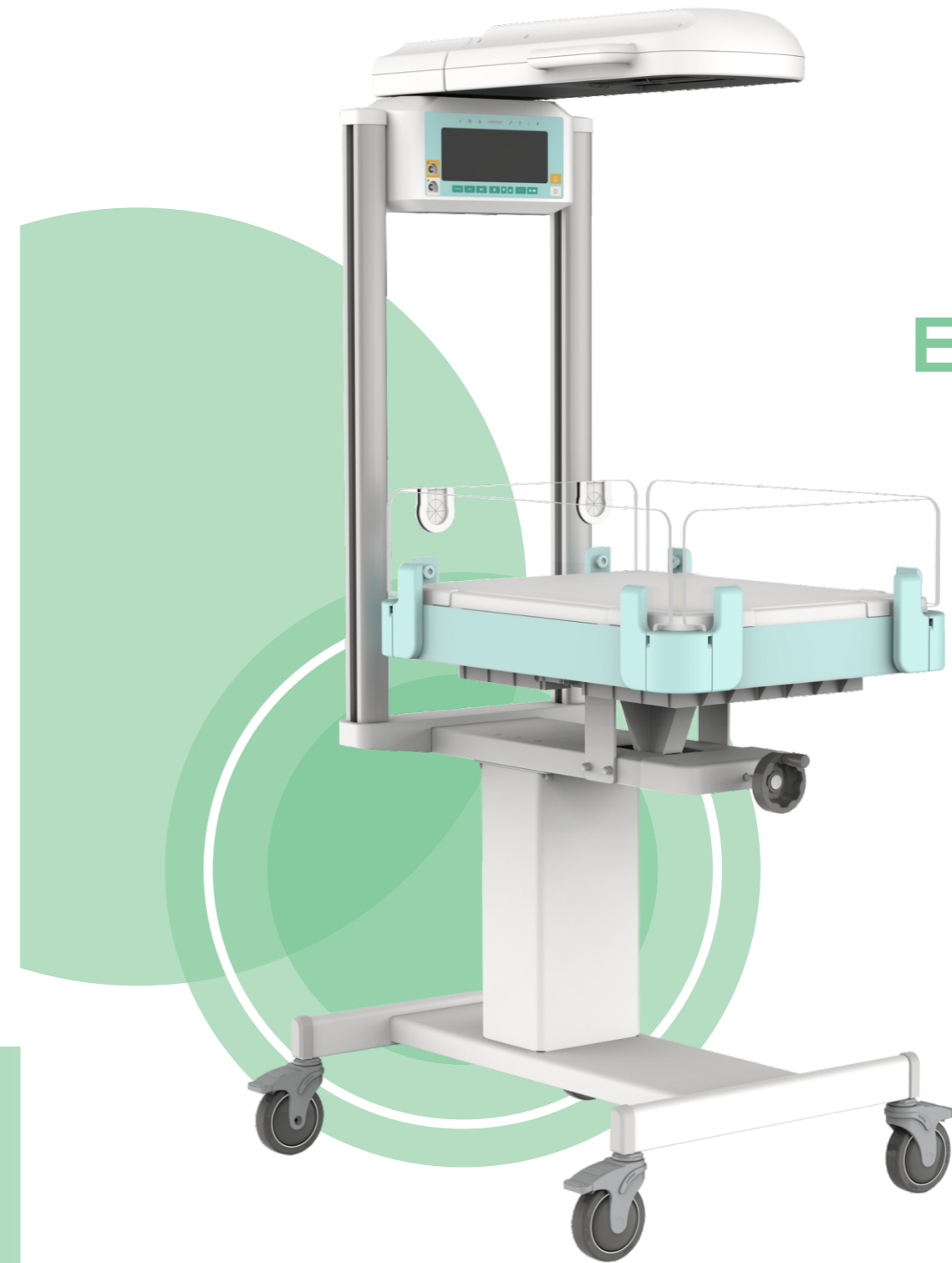

Ecosy630

Infant Warmer



CE ISO
2460 13485

www.bngmed.com

Ecosy 630 INFANT WARMER

Basic configurations

- Far infrared quartz heating tube
- Baby, manual and pre-heating modes
- APGAR and CPR timer
- Transparent guard door on four sides
- Baby bed tilting angle from 0~13°

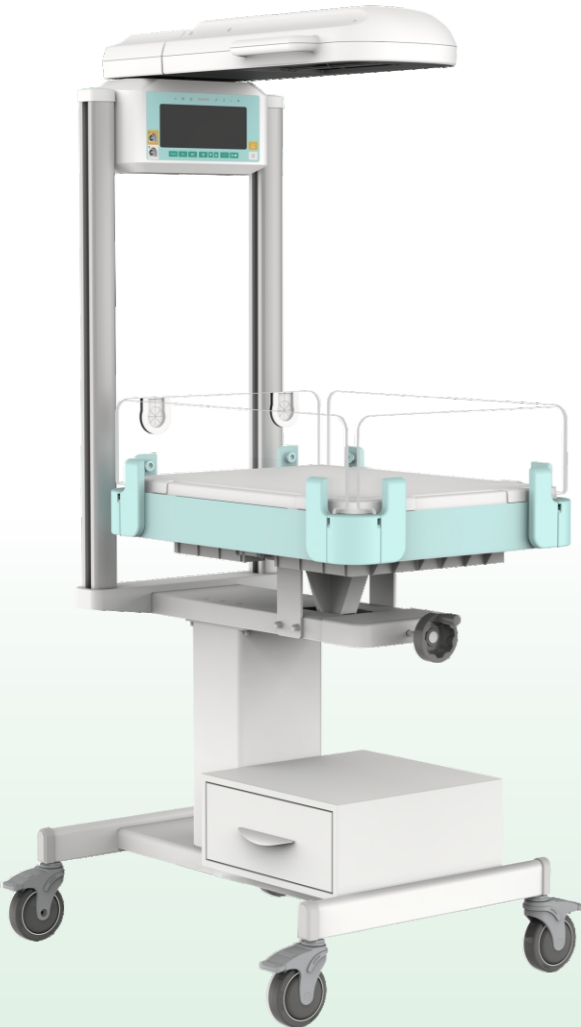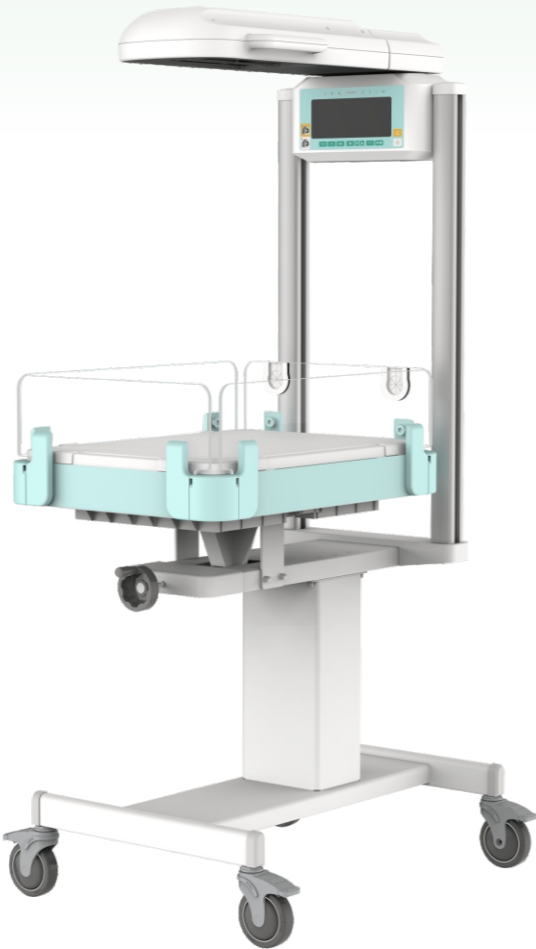


Optional configurations

- Storage drawer
- X-ray cassette
- Lifting stand


■ Heating head module

Heating head adjustable angles: $\pm 120^\circ$ horizontal rotation, $\pm 30^\circ$ vertical rotation.

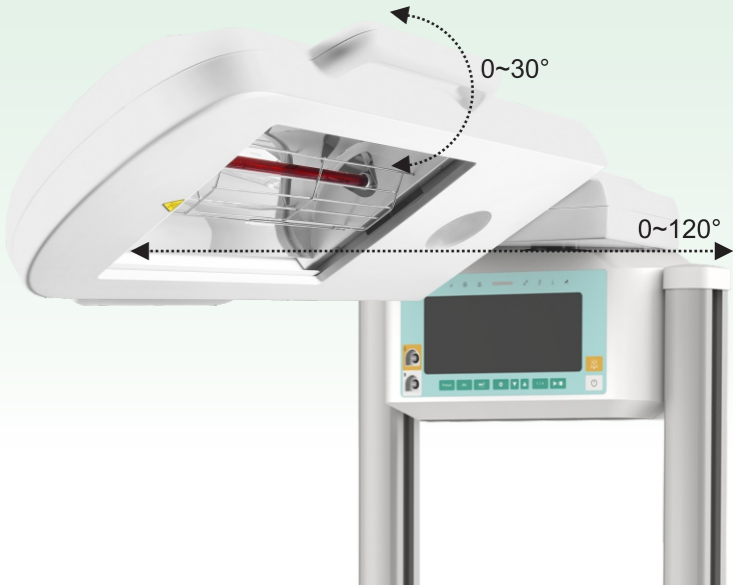


■ LED display

Easily to operate the device and check information.

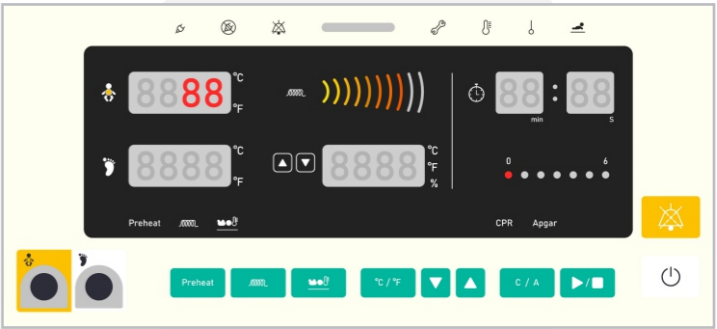



■ Multi-functional Side Rail Support

Design for the external Observation Lamp, LED Phototherapy Unit, Resuscitator Unit, Suction Unit, O2 and Air blender etc.

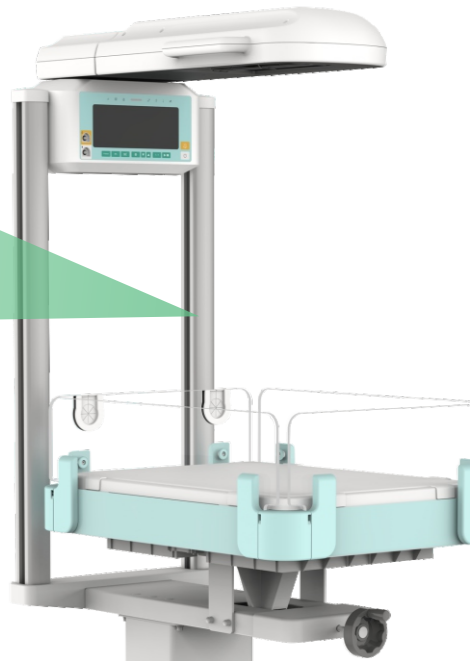

■ **Optional storage drawer**

Storage drawer equipped with damping technology for silent operation and can be operated on both left and right sides which is ideal for storage.

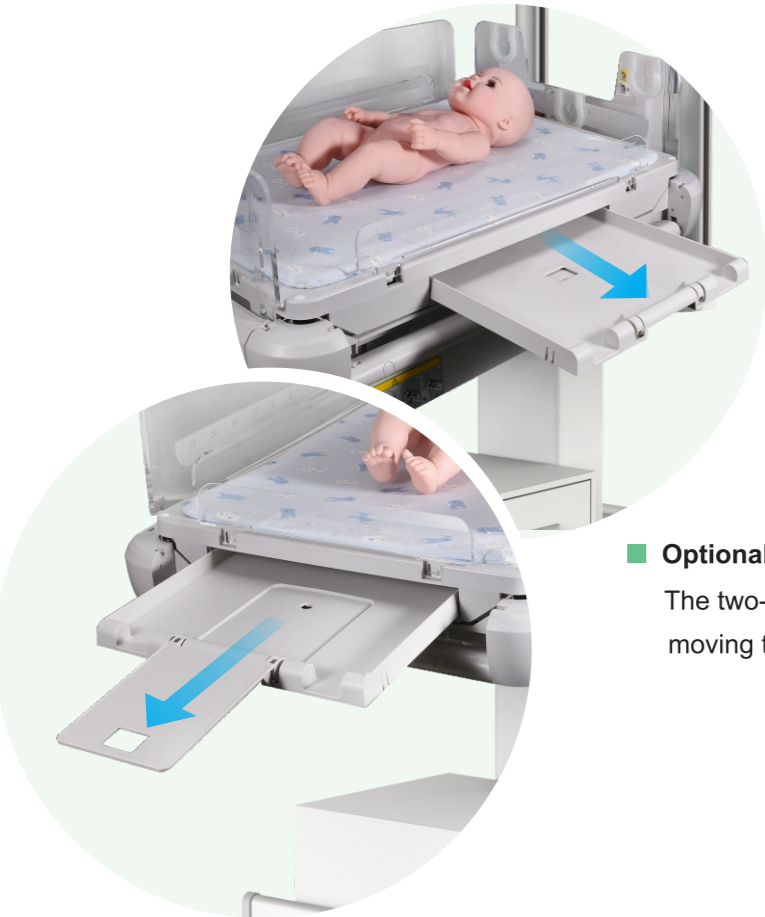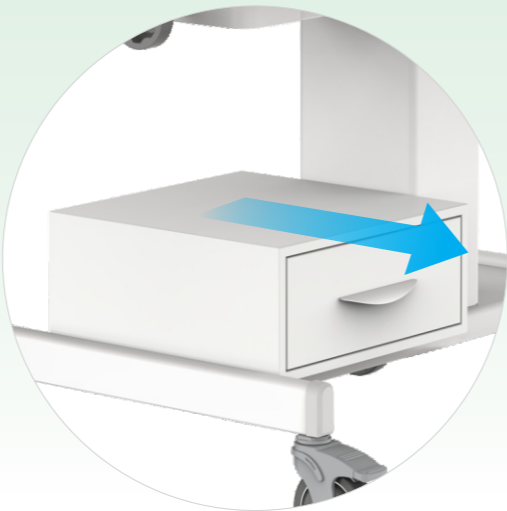

■ **Optional X-Ray cassette**

The two-position pull-out X-ray cassette under the baby tray, avoid moving the newborns when taking the X-Ray examination.

■ **Lifting stand**

Adjust the height of the infant radiant warmer based on operator's height. The foot pedal can be operated from both left and right sides.



Ecosy 630 Infant Warmer

Heater capacity	500W
Temperature control	Servo(Baby) Mode/Manual/Preheating mode
Skin temperature	Setting range:34.0~38.0°C
	Display range:20.0~42.0°C
Heater output	Setting range:0~100%(in 5% increments)
	Indication range:0~100%(in 1% increments)
Timer	APGAR mode/CPR mode switchable
	APGAR Beep:1,5,10 min
Display	LED digital display
Bed tray	Inclination angle(back and forth):±13°
	Mattress surface:500x700x25mm
Dimension	840(L)*1100(W)*1790(H)mm

Optional Functions

X-ray cassette	
Dimension	53x35cm
Height adjustment	Up to 20cm
Storage/Shipping Information	
Temperature	-25~55°C
Humidity	0~93% RH

Notice: Specifications subject to changes without prior notice. All rights reserved by BNG Medical.

